

# SlingBar Mini 220





## Instruction Guide

### Applies to following model:

SlingBar Mini 220  
SlingBar Mini 220 R2R

Prod. no. 3156005  
Prod. no. 3156006

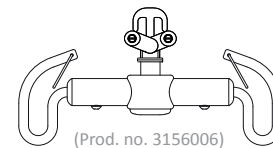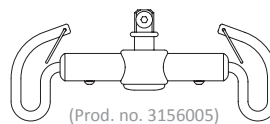

## Product Description

The SlingBar Mini 220 comes in two versions and with the intended use in lifting small patients, for example in a Liko™ TeddySling Low Mod 10. SlingBar Mini 220 (Prod.No 3156005) can be used combined with a Liko™ mobile lift or a Liko™ overhead lift, with fixed assembly. The SlingBar Mini 220 R2R (Prod.No 3156006) version is designed only for Room-to-Room transfers in combination with Liko™ overhead lifts, specified at page 2.

CE SlingBar Mini 220 complies with the requirements of EN ISO 10535.

**Maximum load: 205 kg / 450 lbs.**

### IMPORTANT!

Read the instruction guide for both the patient lift and lifting accessories before use. Lifting and transferring a person always involves a certain level of risk. It is important to completely understand the contents of the instruction guide. The equipment should be used by trained personnel only. Please contact your Hill-Rom distributor in the event of any uncertainties or questions.

## Safety Instructions

Different maximum loads may apply to different products on the assembled lift system: lift, sling bar, sling and any other accessories used. For the assembled lift system, the maximum load is always the lowest maximum load rating for any of the components. Study the markings on the lift and lifting accessories or contact your Hill-Rom representative in the event of questions.



**△ Before lifting, always make sure that the sling's strap loops are correctly fastened to the sling bar hooks when the sling strap is extended, but before the patient is lifted from the underlying surface.**




**△ For safety, load must be applied to all sling bar hooks during the lift!**

## Assembly

**Fixed assembly to:** Liko™ 242 S overhead lift, Golvo™ mobile lift and (Viking™ mobile lift, Uno™ mobile lift)\*

\* An Adapter 12 shall always be installed between the FlexLink and the Slingbar mini 220 to fixed assembly of Viking™ mobile lift and Uno™ mobile lift, see below

Tools required: 1 x  13 mm (1/2 inch.)  
1 x  8 mm (5/16 inch.)

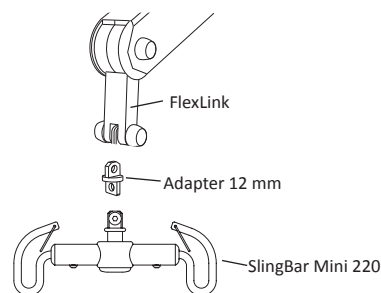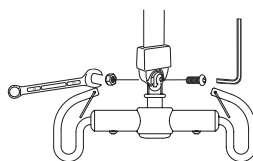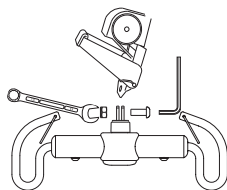

**Adapter 12 mm,**  
See applicable assembly instruction!


Prod. No. 2016504



**Fixed assembly to:** LikoGuard™ overhead lift

Require, Fixed kit 4,  
**NOTE! Shall be assembled by personnel authorized by Hill-Rom!**

Prod. no. 3308860

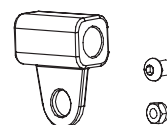


Continues on next page 


### Attaching the SlingBar Mini 220 R2R to:

Likorall™ 242 R2R overhead lift, Likorall™ 243 overhead lift, Likorall™ 250 overhead lift, Likorall™ 200 overhead lift.

The SlingBar Mini 220 R2R double hook with its Quick-release function, has protection against unintentional detachment, see fig. 1.

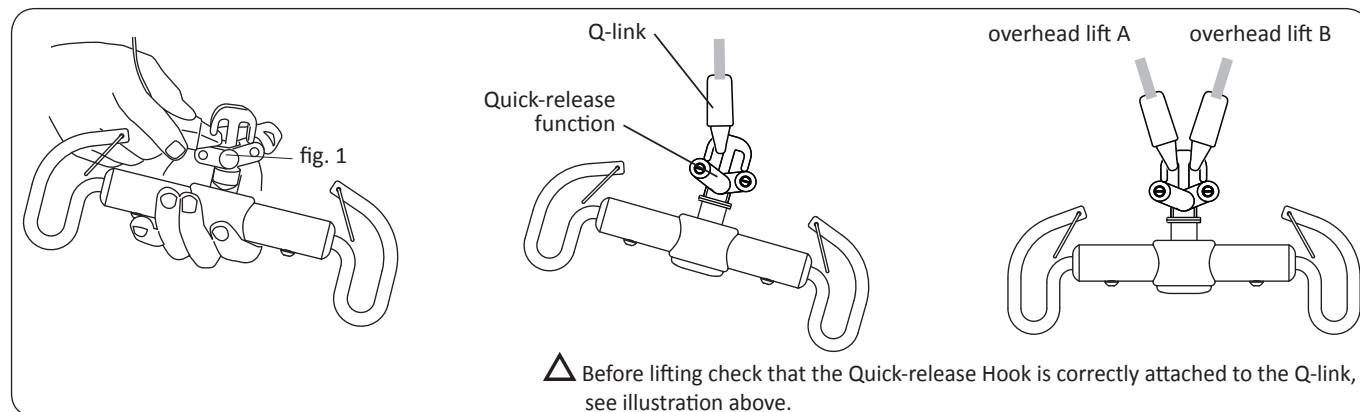

### Installation of latches

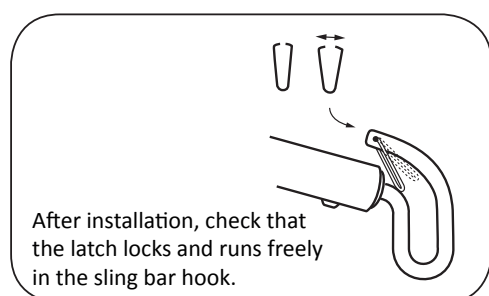

## Care and maintenance

For safe and trouble-free operation, a few routine procedures should be performed every day, before the sling bar is in use:

- Check the screw and locking nut that attaches the sling bar
- Check the function of latches and that the latches are intact; missing or damaged latches must always be replaced
- When using a Quick-release Hook system, check that the Quick-release Hook is correctly fastened to the lift and the sling bar
- Before lifting, check that the lifting accessory hangs vertically and can move freely.

When necessary, clean the lift with a moist cloth. Find more detailed information regarding cleaning and disinfection of your Liko™ product in the document *Care and Maintenance* at our website: [www.liko.com](http://www.liko.com).

### Service

A periodic inspection of the sling bar should be carried out at least once a year.

△ **Periodic inspection, repair and maintenance may be performed only in accordance with the Liko™ service manual by personnel authorized by Hill-Rom and using original Liko™ spare parts.**

### Expected Life Time

The product has an expected life time of 10 years when correctly handled, serviced and periodically inspected in accordance with Liko's instructions.

### Recycling

For instructions on how to recycle your sling bar, please visit our website: [www.liko.com](http://www.liko.com).

### Product Changes

Liko's products undergo continuous development, which is why we reserve the right to make product changes without prior notice. Contact your Liko representative for advice and information about product upgrades.

### Design and Quality by Liko in Sweden

Hill-Rom's Management system are certified in accordance with ISO 9001 and its equivalent for the medical device industry, ISO 13485.

Hill-Rom's Management system is also certified in accordance with environmental standard ISO 14001.



[www.hill-rom.com](http://www.hill-rom.com) | [www.liko.com](http://www.liko.com)