


Hill-Rom® OB-GYN Stretcher

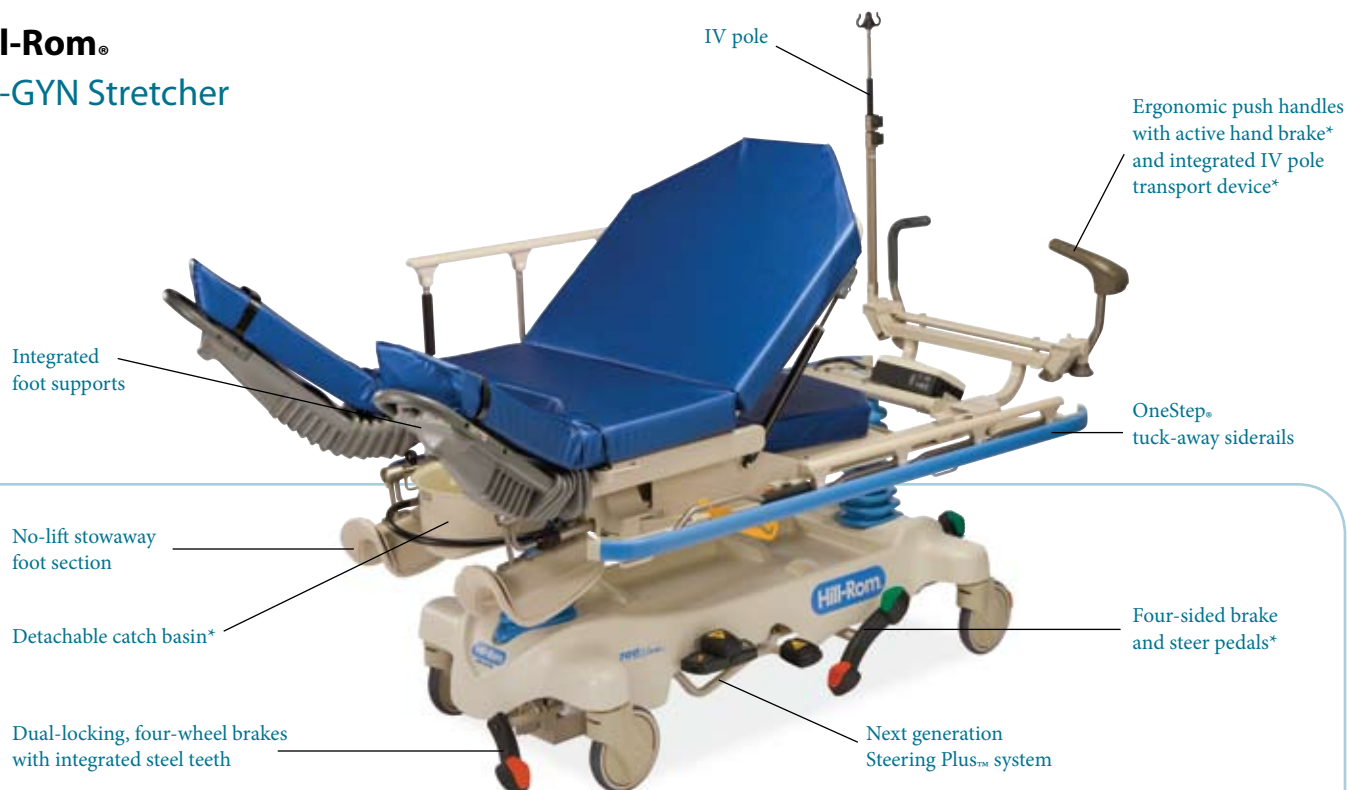

At a glance

Hill-Rom stretcher key features:

- **700 lb weight capacity**
- **Sliding litter** enables quick and easy patient positioning
- **Stowaway foot section** eliminates heavy lifting and storage concerns
- **Integrated foot supports** save time by eliminating the need to search for stirrups and calf supports
- **Integrated catch basin***, combined with patient grip handles, offers convenience and saves time
- **Next generation Steering Plus™** system for crisp cornering and control
- **Dual-locking four-wheel brakes** for a strong, stable platform
- **Power washable**
- **Three year parts warranty** and lifetime frame weld warranty

* Optional feature

Hill-Rom® OB-GYN Stretcher

With its fast and easy patient positioning, Hill-Rom's OB-GYN stretcher (P8050) is ideal for use in both maternity and emergency. Built to last year after year, this stretcher offers outstanding mobility, convenience and quality.

Mobility. Why struggle to steer and stop? An optional active hand brake reduces strain of going down ramps and sudden stops. Our integrated IV pole transport device helps you to maneuver stretcher and IV pole stands with ease. The next generation Steering Plus™ system allows crisp cornering and control.

Convenience. Save time with this stretcher's integrated foot supports, sliding litter and optional catch basin. Reduce heavy lifting and storage concerns with our stowaway foot section. An optional built-in fiber optic exam light eliminates searching for another piece of equipment.

Quality. Our reinforced frames are made of heavy-gauge, automotive grade steel for long-lasting quality — and a 700 lb weight capacity to handle higher patient weights. Nickel-plated push handles, IV poles, siderails and other components offer superior rust resistance, even in harsh environments.

Hill-Rom®

Enhancing Outcomes for Patients and Their Caregivers.™

Hill-Rom®

OB-GYN Stretcher

Features	Standard Features	Optional Features
700 lbs weight capacity	•	
Sliding patient litter	•	
Integrated foot supports	•	
No-lift stowaway foot section	•	
Active hand brake		•
Integrated IV pole transport device		•
Brake/Steer pedals – head and foot	•	
Brake/Steer pedals – all four sides		•
Next generation Steering Plus™ system	•	
Ergonomic push handles		•
Permanent IV pole – 3 stage (head only)		•
30-inch width	•	
Siderail system	Tuck-away	
Radiolucent deck	•	
3-inch high density foam mattress with fire barrier	•	
Factory mounted light		•
Bumper colors – choice of neutral, blue, red, teal or purple	Choice	



Integrated Foot Supports



Integrated IV Pole Transport Device



Active Hand Brake



Factory Mounted Light

Technical Specifications

Maximum weight limit.....	700 lbs (317.5 kg)
Low position – floor to top of surface.....	24.5" (62.2 cm)
High position – floor to top of surface.....	37" (94 cm)
Overall length	80.5" (204.5 cm)
Overall width (siderails up)	35.5" (90.2 cm)
Overall width (siderails stored)	31" (78.7 cm)
Siderail length.....	47" (119.4 cm)
Siderail height above mattress/surface.....	13" (33.0 cm)
Mattress size	30" (76.2 cm) x 75" (190.5 cm)
Maximum head elevation	90°
Maximum Trend/Reverse Trend	18°
Floor to base clearance	3.5" (8.9 cm)
Caster size.....	8" (20.3 cm)

Accessories

- Stretcher pillow
- IV transporter
- Patient tray
- Utility tray
- Oxygen tank holder
- Liquid oxygen tank holder
- Infusion support system
- Siderail covers
- Armboard
- Detachable headboard
- Lateral x-ray cassette holder
- Chart holder
- Convertible footboard/headboard
- Accessory IV pole
- Transport straps
- Foot extender
- Paper roll dispenser
- Integrated telescoping calf supports
- Integrated catch basin

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

© 2009 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.

145650 rev 2 8/26/09



Enhancing Outcomes for Patients and Their Caregivers.™

USA 800-445-3730
Canada 800-267-2337

www.hill-rom.com